


Art nr: 15471, 15472

Ingående delar

1 x

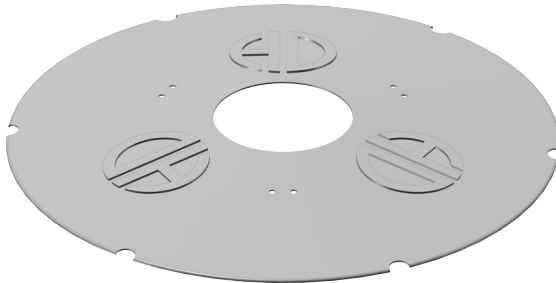

3 x




Montering


1



2

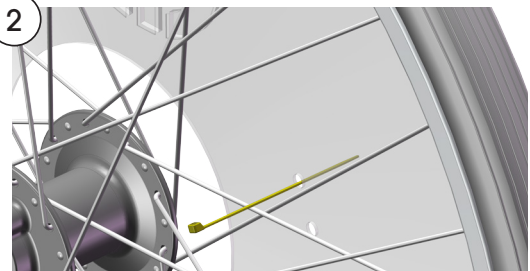

3

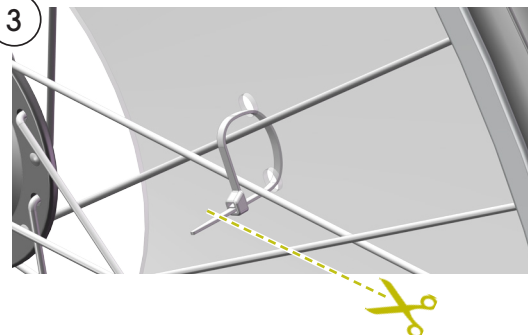

Ta bort skyddsplasten från båda sidor av ekerskyddet.

1. Passa in ekerskyddet mot ekrarna på hjulets utsida genom att trä det innanför drivringarna. Se till att hålbilden stämmer med ekermönstret och alla hålpar träffar ekrar. Roter annars till rätt position.
2. Montera de medföljande buntbanden från hjulets insida enligt figur. Buntbandens låshuvud ska hamna på insidan.
3. Dra åt buntbanden och klipp av de utstickande ändarna.

Skötselråd

Ekerskydd kan rengöras genom avtorkning med mildt rengöringsmedel och vatten. Byt alltid ut skadade delar omedelbart. För övrig information se Underhållsanvisning, art. nr. 95730.

Spoke protection


User manual / Assembly instructions

Page 1 (1)

Art. no.: 15471, 15472

Parts included



1 x

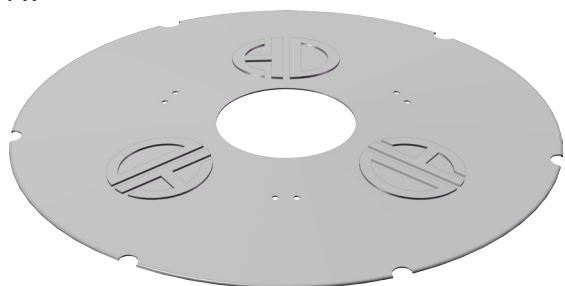

3 x



Assembly

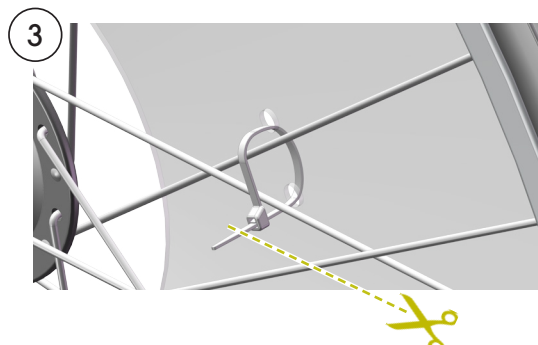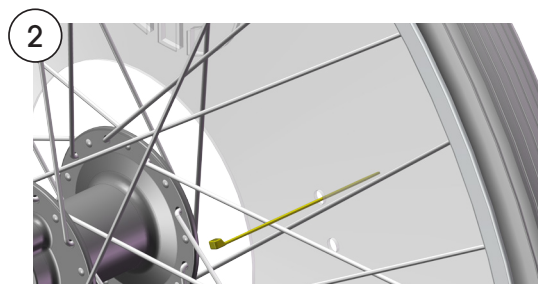

Remove the protective plastic film from both sides of the spoke protector.

1. Fit the spoke protector against the spokes on the outside of the wheel behind the hand rim. Make sure that the holes align properly with the spokes, otherwise rotate the protector to correct the positioning.
2. Fit the cable ties from the inside of the wheel as shown. The locking head of the cable should be on the inside.
3. Pull the cable ties tight and cut off the extra material.

Care instructions

The spoke protector can be cleaned by wiping with mild cleaning agent and water. Always replace damaged parts immediately. For more information see the Maintenance Instructions, art. no. 95730-1.